Antibacterial **Bioelectric Bandage Kit Instructions**

PowerHeal™

EN

PowerHeal™ Antibacterial Bandages feature embedded microcell batteries that generate a low level of electricity when in the presence of a conductive medium such as water-based wound hydrogels, saline, or water. Note: Products are not made with natural rubber latex.

PowerHeal™ Bioelectric Bandage with Adhesive has 3 layers: a broad-spectrum antibacterial wound contact layer, an absorbent lyocell layer, and a polyurethane adhesive layer to hold the bandage in place and maintain a moist wound environment.

PowerHeal™ Bioelectric Bandage with No Adhesive is a broad-spectrum antibacterial polyester single wound contact layer. It requires an alternate covering and optional absorbent pad (sold separately) to secure it and maintain a moist wound environment.


Take charge of your healing.

PowerHeal.com


How To Use

Bioelectric Bandage with Adhesive


Please Note: CLEANSE wound with an appropriate wound cleanser. Allow area to dry completely before applying the bandage. This will help prevent skin irritation and ensure the adhesive sticks optimally.

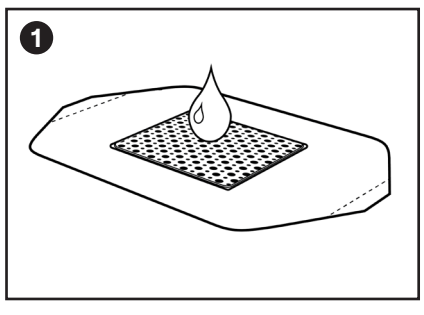Moisten the dotted bandage with a **thin, even layer of Energel® wound hydrogel.**

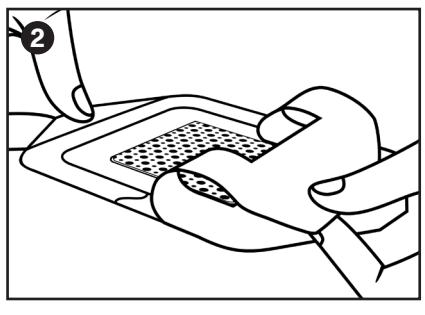
Peel off the paper liner and apply the dotted surface to your wound. Gently press adhesive down to firmly seal.

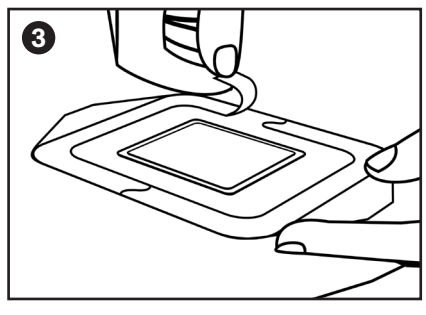
Remove the paper frame on the top side of the bandage and gently press to secure the adhesive border.

Bioelectric Bandage with No Adhesive


Please Note: CLEANSE wound with an appropriate wound cleanser. Allow area to dry completely before applying the bandage. This will help prevent skin irritation and ensure the adhesive sticks optimally. Prior to moistening, PowerHeal™ may be cut to fit if needed. Cut it large enough so that it will extend ½ - 1 inch beyond the wound edges (cutting may result in slight fraying of the bandage edges).

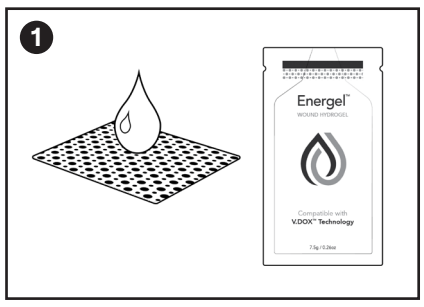Moisten the dotted side of the bandage with a **thin, even layer of Energel® wound hydrogel.**

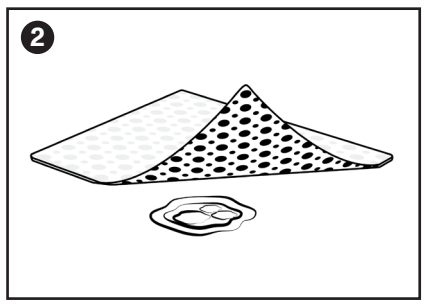Apply the dotted side down on the wound. Gently press to ensure direct contact with the wound.

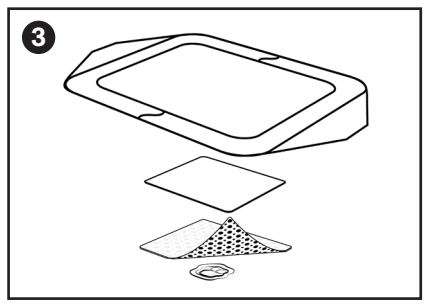Cover the bandage with your preferred sterile secondary dressing (i.e. absorbent pad, gauze, medical wraps) to secure in place and maintain a moist wound environment.

Indications

PowerHeal™ Antibacterial Bandage for over-the-counter use is intended for the management of wounds, to provide a moist wound environment, and is indicated for superficial wounds such as minor cuts, scrapes, irritations, abrasions, blisters, etc.

Energel® Wound Hydrogel is used to cover a wound to protect against abrasion, friction, drying out, and contamination. It can also be used to activate PowerHeal™ products.

Contraindications

PowerHeal™ Antibacterial Bandages—Do not use on individuals with sensitivity or allergy to silver, zinc, or other bandage components. (Elemental silver and elemental zinc are battery components used in the bandage to preserve it and minimize or prevent the growth of bacteria within the bandage.)

Energel® Wound Hydrogel—None known.

REF	Catalog number
LOT	Batch code
	Use by date
	Manufacturer
STERILE R	Sterilized using irradiation
STERILE EO	Sterilized using ethylene oxide
LATEX	Product is not made with natural rubber latex
	Do not use if package is damaged and consult instructions for use
2	Do not re-use
MR	(MR) Unsafe
	Consult instructions for use or consult electronic instructions for use

Precautions—PowerHeal™ Bandages

- Do not use PowerHeal™ Bandages with topical agents such as antibiotic creams or ointments, petroleum-based products, enzymatic debriders, silver or zinc creams, or oxidizing agents.
- Remove PowerHeal™ during energy-based procedures (such as radiofrequency, ultrasound, or radiation) where the bandage may interfere with delivery.
- MR unsafe. Remove PowerHeal™ prior to a MRI or HBOT procedure and apply a new bandage after the procedure.
- Avoid direct bandage contact with electrodes or conductive gels during electronic measurements; e.g., EEG (electroencephalogram) or ECG (electrocardiogram).

Duration of use:
- Single use only. For external use only.
- Do not use if the pouch is open or damaged prior to use.
- The safety of daily PowerHeal™ Bandages use for longer than 28 days has not been studied.

Uncontrolled bleeding, infection, sensitive skin, allergies:
- Stop using PowerHeal™ and consult a physician if allergy, irritation, increased pain, maceration, or irregular skin discoloration occurs.
- PowerHeal™ is not intended to be used on wounds with uncontrolled bleeding.
- Infected wounds should be treated under a physician's supervision.
- The use of adhesive bandages on fragile or sensitive skin may pose a risk of skin damage upon removal of the bandage. To help minimize this risk, a skin prep may be applied on the wound perimeter prior to bandage application (not included).

Warning—PowerHeal™ Bandages

Frequent use of this product may result, in rare occasions, temporary discoloration of the skin.

Precautions—Energel® Wound Hydrogel

- Single use only. For external use only. Discard sachet after single use.
- Sterile until opened.
- Discard product if packaging is damaged or evidence of contamination is present.
- Bleeding should be controlled before application of a bandage.
- Do not use if you have a sensitivity to parabens.

How To Take Care of Your Bandage

Bandages should be changed based on the condition of your wound and the surrounding skin. For example, if you have a lot of drainage coming from your wound, you may need to change it once or twice a day, especially in the early stages of healing. As your wound progresses in its healing, it is typical to see less drainage. In this case, you may want to change your bandage every 2 to 3 days. PowerHeal™ bandages may be left in place for up to 7 days.

You may shower with the PowerHeal™ bandage in place. Avoid direct spray of water onto the bandage to avoid saturation and premature adhesive release.

How To Remove Your Bandage

To remove your PowerHeal™ bandage, lift one corner and slowly pull the bandage back, keeping it parallel with the skin surface until it is fully removed. **Please Note: If the bandage sticks to your wound, do not force it off. Instead, moisten the bandage and wound area thoroughly with clean water or saline and let it soak in until the bandage can be removed without tissue disruption.**

PowerHeal Wounds, LLC
Tempe, AZ 85288, USA
(866) 496-8743 | www.powerhealwounds.com

© 2023 Vomaris. All rights reserved.
PowerHeal is a trademark and Energel is a registered trademark of Vomaris Innovations, Inc. PATENTED. See www.vomaris.com.

IFU-118 Rev A

Prefer to watch a demo instead?
Scan this QR code to view our video demo.

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